

आनंद रंगपिल्लैग्रंथालय / ANANDA RANGAPILLAI LIBRARY

पांडिच्चेरी विश्वविद्यालय / PONDICHERRY UNIVERSITY

Kalapet, Puducherry - 605014

Dr.M Vijayakumar
University Librarian

☎ : + 91 - 413 2654205
librarian@pondiuni.ac.in

Ref: No. PU/ARPL/JLS/2026/

23rd September 2026

E-CIRCULAR

Trial Access to SYSTAT 13.2 (Statistical Analysis and Graphics Software)

The Ananda Rangapillai Library, Pondicherry University, is pleased to inform you that, as a continuation of the Hands-on Training on *SYSTAT13.2 – Powerful Statistical Analysis and Graphics Software* conducted on 03 September 2026, a 30-day trial access to *SYSTAT 13.2* has been made available for the benefit of the Pondicherry University community.

SYSTAT 13.2 is a powerful statistical computing platform that supports statistical analysis, data management, hypothesis testing, statistical modelling, advanced multivariate techniques, and scientific data visualisation. Faculty members, research scholars, and students are encouraged to utilise the trial access of the *SYSTAT 13.2* for their academic, teaching, research, and data analysis requirements.

The *SYSTAT 13.2* software can be downloaded using the link given below:

<https://inpixon.files.com/f/8b6d93dc0b01746a>

Serial Number: 797791234

The *SYSTAT 13.2* installation guide is also attached to this e-circular. Users are requested to follow the installation instructions carefully while installing and activating the software.

Important instructions:

- The trial access is available for 30 days.
- The software is intended strictly for use by the Pondicherry University user community.
- Users are requested to install the software and complete the activation process by following the attached installation instructions.
- The serial number provided above should be used only for the authorized trial access.

For more details, please contact the Journals Section of Ananda Rangapillai Library, Pondicherry University, Email: journal.lib@pondiuni.ac.in, Phone: +91-413 2654282.

Vijayakumar


Dr. M. Vijayakumar
University Librarian

To:

- The Head, Computer Centre, PU – with a request to host this circular on the Pondicherry University website.
- The In-Charge, DLS, Ananda Rangapillai Library – with a request to host this circular on the Library website.


SYSTAT® 13 Installation Instructions

Installing SYSTAT

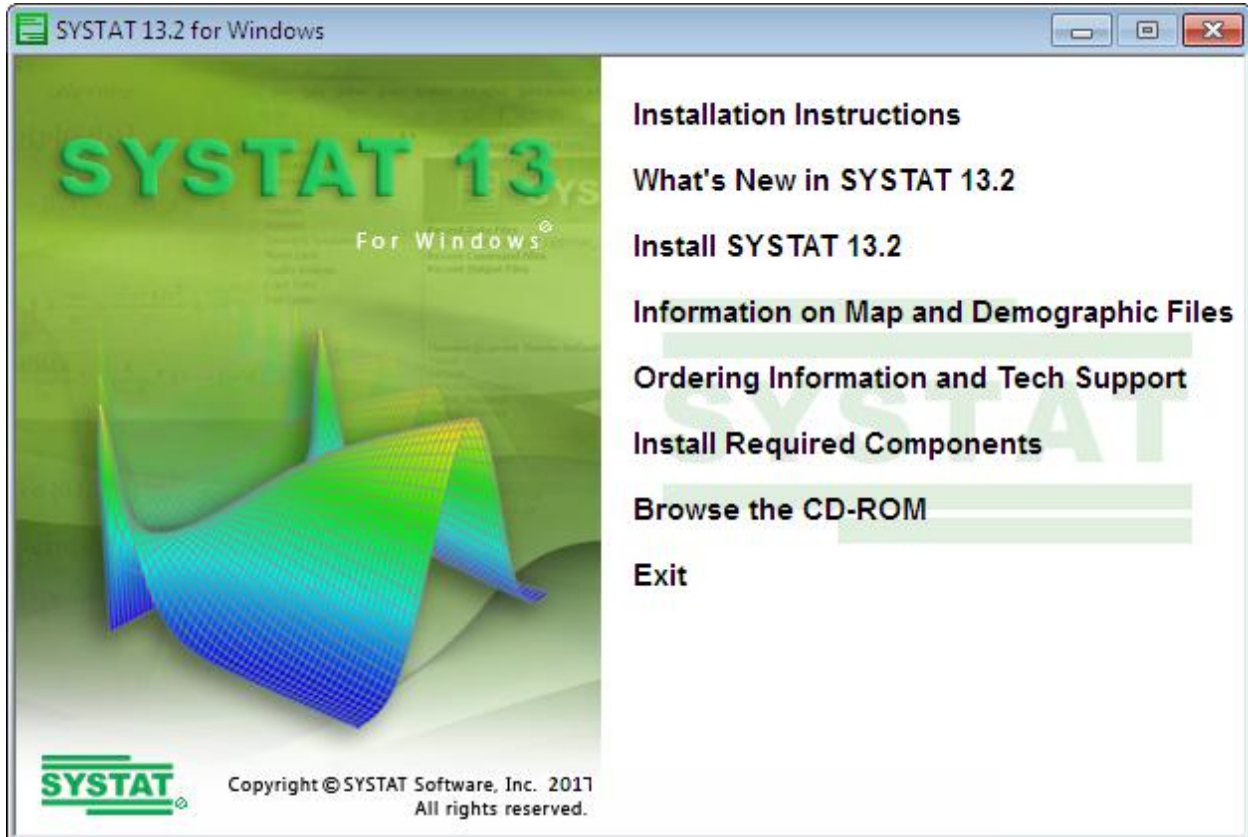

To install SYSTAT, insert the SYSTAT CD into the CD-ROM drive. The AutoPlay feature will present a special menu. To run the Setup program, choose **Install SYSTAT 13**, and then follow the instructions that appear on the screen. In case the Autorun is disabled then open the CD from Windows Explorer and click on the Autoplay.exe file.

The Setup program installs SYSTAT 13 into a location of your preference specified during the installation. At the beginning of the installation, you are prompted to enter your name, the name of your company and the serial number (printed on the back of the product CD's jewel case). This information is later displayed in the SYSTAT 13 Help/About box for your reference. Please note that filling the serial number field is mandatory for the installation to go forward.

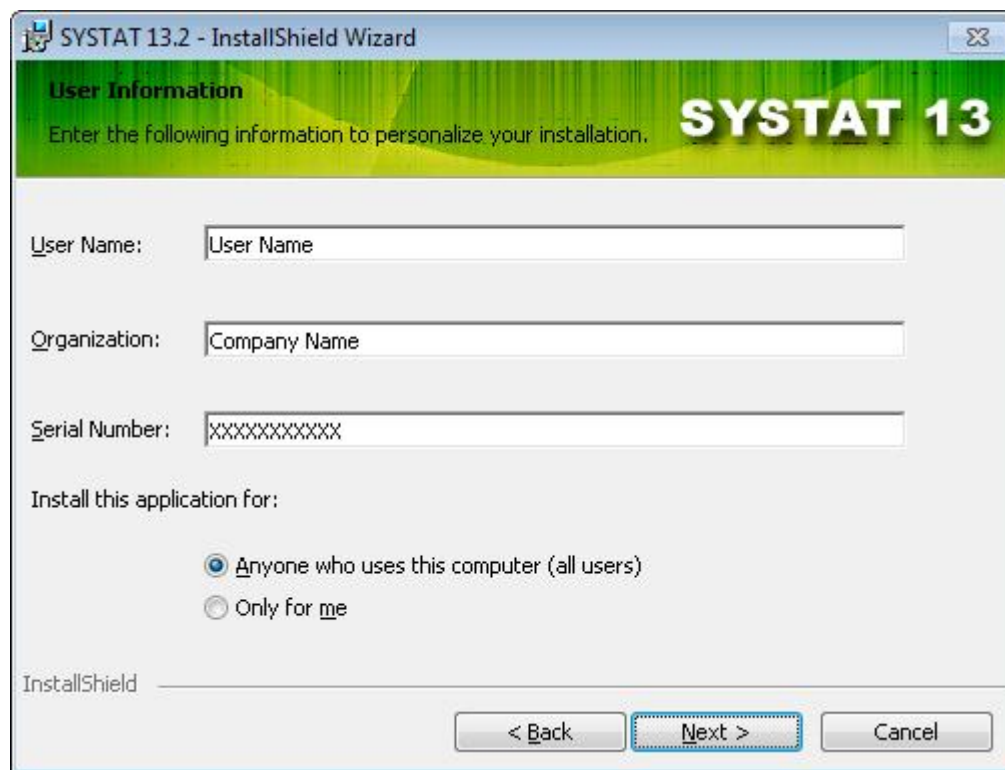

The image shows a Windows-style installation wizard window titled "SYSTAT 13.2 - InstallShield Wizard". The window has a green header bar with the text "User Information" and "Enter the following information to personalize your installation." followed by the "SYSTAT 13" logo. Below the header, there are three text input fields: "User Name:" with the placeholder "User Name", "Organization:" with the placeholder "Company Name", and "Serial Number:" with the placeholder "XXXXXXXXXX". Below these fields, the text "Install this application for:" is followed by two radio button options: "Anyone who uses this computer (all users)" (which is selected) and "Only for me". At the bottom left, the text "InstallShield" is visible. At the bottom right, there are three buttons: "< Back", "Next >" (which is highlighted with a blue border), and "Cancel".

The above installation creates the following shortcuts in the Start menu:

Start => Programs => SYSTAT 13 =>

Command Translator

Manuals

Restore Default Settings

SYSTAT 13

SYSTAT FAQ

SYSTAT Help


What's New in SYSTAT 13

Uninstall

Using the next screen a common destination folder for all features can be selected



The above dialog box lists the available features in the installation. To prevent any of the features from being installed click on that feature and select "Entire feature will be unavailable." item as shown in the following dialog box

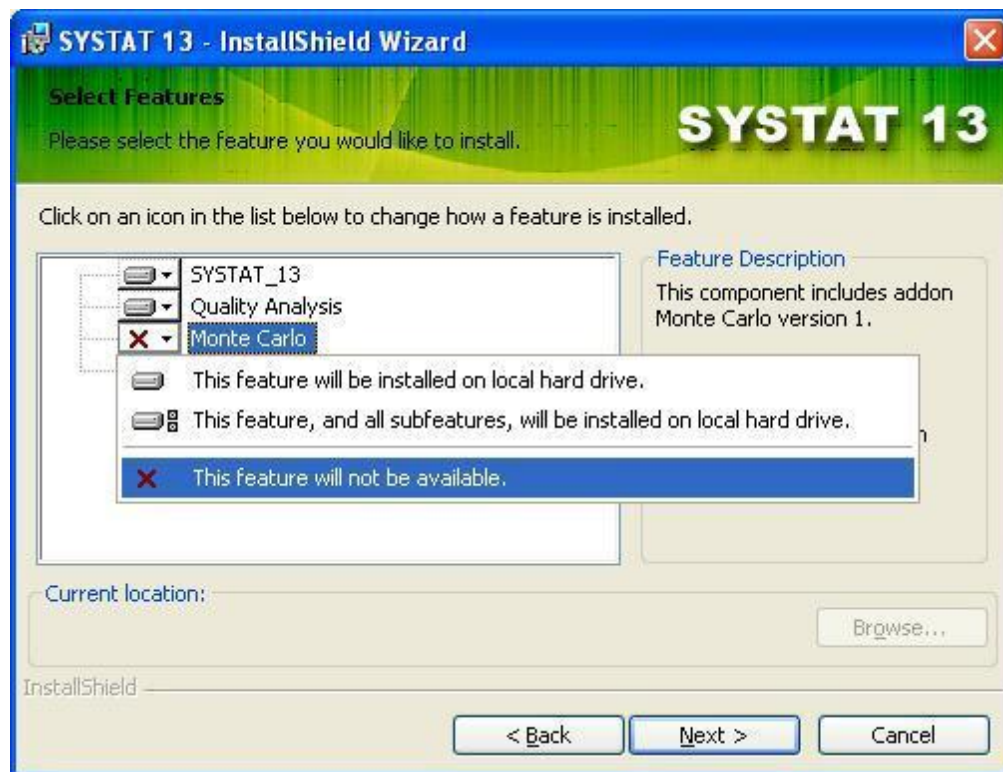

At the end of installation the following dialog box appears; by default SYSTAT 13 manuals are slated to be installed; to prevent installation of SYSTAT 13 manuals uncheck the button 'Install SYSTAT pdf Manuals'.



Obtaining a SYSTAT License

The standalone and network installations contain 30 days evaluation license. The license manager for network installation also contains 30 days evaluation license for 2 users.

Updating SYSTAT License

Systat Application startup will prompt for License update or Application menu "Help -> Update License". The following dialog box opens.

Systat 13.2 License Utility

Help

License Activation Log Info License Information Contacts

Available License

Product Name	Version	License Type	Expiry in Day(s)	License Manager
Systat (ExactTest)	13.2	Trial	30	WIN7VM38-64

Thank you for evaluating Systat. You can upgrade this to a full version at anytime by contacting one of the sales offices listed in the information tab. If you already have a new license, pick the new license file and click the 'Apply License' button below to install it.

New License

☒ Product Key ☐ V2C File

Product Key

License File Path Browse....


Apply License

Close

If you have a product key, Select Product Key radio button and enter value in the Product Key edit box.

If you have a V2C license file then Select V2C File radio button then browse to your license file or you can type the absolute path to your license file in the license file edit box.

Next click the Apply License button. When the Apply License button is clicked License gets updated.

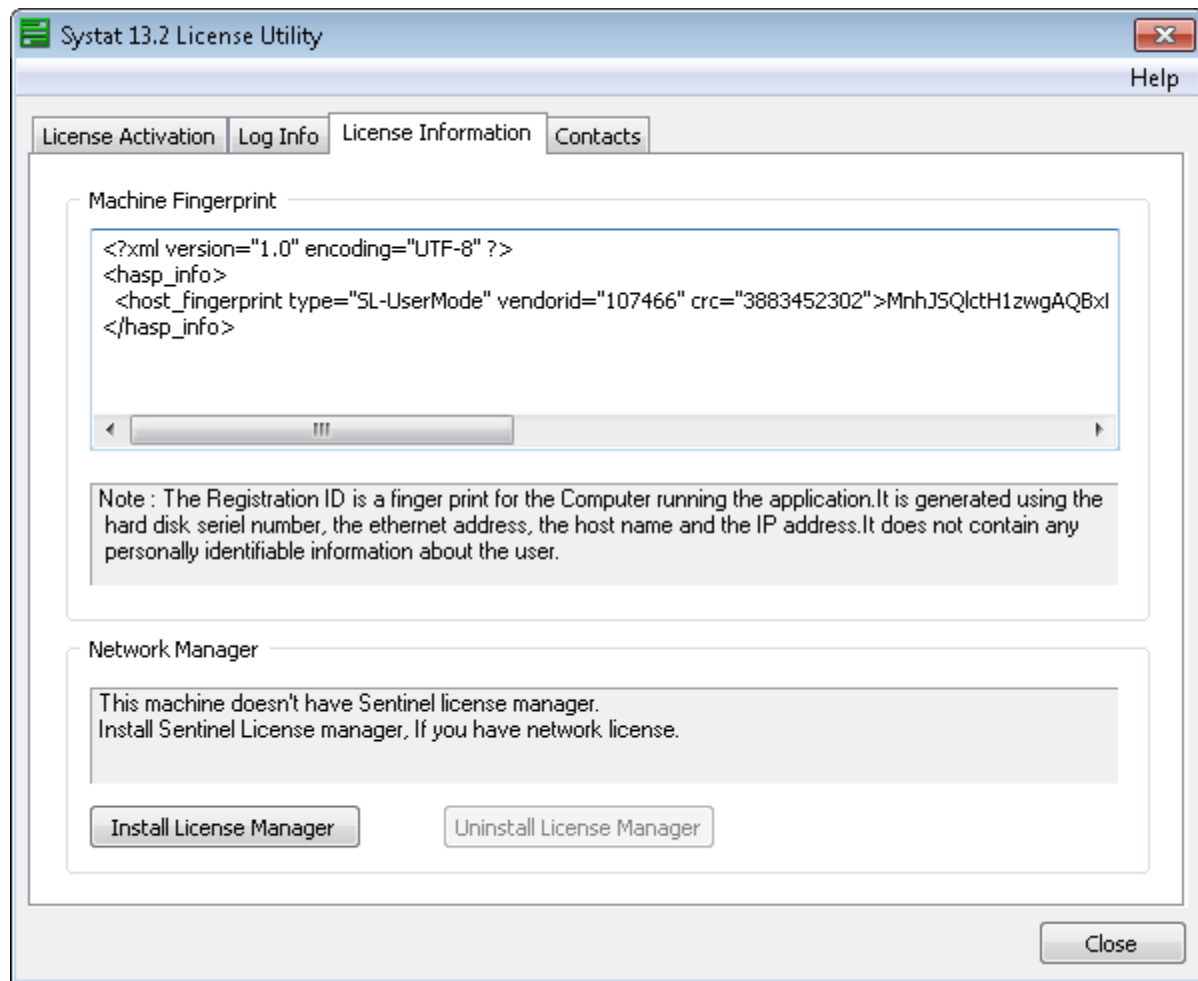


Silent installations

Systat13 provides an option to install the application silently using system `msiexec.exe`. There are three ways, Systat13 can be installed silently.

- Using `SystatSilentInstallation.exe` by reading `SilentInstallProperty.ini` property file.
- Using system's `msiexec.exe` (at command prompt) by reading `SilentInstallProperty.ini` property file.
- Using system's `msiexec.exe` by directly providing properties command prompt

About SilentInstallProperty.ini file

Property file specifies the parameters to be passed to systat13ENU32.msi to have end user select options that are available at User Interface level. Please refer below the SilentInstallProperty.ini content and its description.

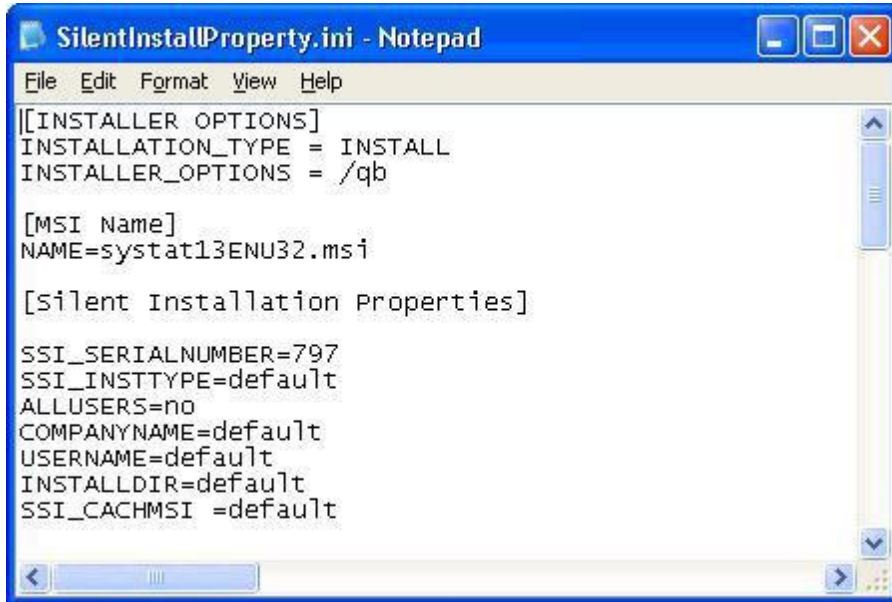

```
[[INSTALLER OPTIONS]
INSTALLATION_TYPE = INSTALL
INSTALLER_OPTIONS = /qb

[MSI Name]
NAME=systat13ENU32.msi

[Silent Installation Properties]

SSI_SERIALNUMBER=797
SSI_INSTTYPE=default
ALLUSERS=no
COMPANYNAME=default
USERNAME=default
INSTALLDIR=default
SSI_CACHMSI =default
```

[Details Description about each property]

Installation Properties	
SSI_SERIALNUMBER	Represents a serial number of systat13 application. This is mandatory property. if you do not provide serial numbers, application will not be installed.
SSI_INSTTYPE	Represents type of license to be used for systat 13 application. This is optional property. *default=Standalone other valid values= Network or Site
ALLUSERS	Represents type of application installation i.e. Per Machine Installation or Per User Installation. This optional property. ALLUSERS=yes(1), to install application for all users i.e. Per Machine Installation. Not Recommended for Silent Installation. Note: This option to work, user should have administrative

	privileges to install the application, otherwise it will throw an error. ALLUSERS=no(""), to install the application for current user. This is recommended option during silent installation. *default= yes
INSTALLDIR	Represents application's target directory to which application files will be copied. *default= "SYSTEM DRIVE:\Program Files\SYSTAT 13" i.e "C:\Program Files\SYSTAT 13" Other valid values= "Any valid path" such as "D:\Softwares\SYSTAT"
COMPANYNAME	Represents name of the organization. *Default = Organization name registered to Operating System. Other Valid Value= "Enter your name of Company"
USERNAME	Represents name of the user to whom operating System has been registered *default=User Name registered to Operating System. Other Valid Value= "Enter Your User Name"
SSI_CACHMSI	Represents cache MSI locally for future maintenance. *default=yes, to cache the MSI and .cab files locally for future maintenance SSI_CACHMSI=no, do not cache the MSI and .cab files locally.
Below options are valid only when Systat 13 is installed silently using SystatSilentInstallation.exe	
INSTALLATION_TYPE	INSTALL - To install the Systat 13 application UNINSTALL – To uninstall the already installed Systat 13 application.
INSTALLER_OPTIONS	It can be either of following /qb - displays basic UI with progress /qn - Displays No UI /qr - Reduced UI

[MSI Name]	Provide the name of the msi to be launched. If you are installing the 32bit application, the name would be systat13ENU32.msi and for 64bit, the name would be systat13ENU64.msi *default - if you provide Property value as 'default', msi would initialize values defined in Installer Property Manager.
------------	---

Using SystatSilentInstallation.exe by reading SilentInstallProperty.ini property file.

Copy the Systat13ENU32.msi (along with SYSTAT~1.cab, Qualit~1.cab, Monte_~1.cab and ExactT~1.cab), SilentInstallProperty.ini and SystatSilentInstallation.exe from CD ROM to C drive. Modify the SilentInstallProperty.ini according to requirements and double click on the SystatSilentInstallation.exe. Depending on INSTALLER_OPTIONS (/qb or /qn) in SilentInstallProperty.ini, either basic UI with progress bar or with No UI will be displayed.

Using system's msixec.exe (at command prompt) by reading SilentInstallProperty.ini property file.

For installing Systat13 through command prompt silently by reading the SilentInstallProperty.ini property file, copy the Systat13ENU32.msi (along with SYSTAT~1.cab, Qualit~1.cab, Monte_~1.cab and ExactT~1.cab) and SilentInstallProperty.ini from CD ROM to C drive, open the command prompt and run the below command. Please make sure that SilentInstallProperty.ini has been changed according to needs before executing the command.

```
msiexec /i C:\Systat13ENU32.msi /qb (for 32bit Application)
msiexec /i C:\Systat13ENU64.msi /qb (for 64bit application)
```

Using system's msixec.exe by directly providing properties command prompt

For installing Systat13 through command prompt silently by providing properties at command prompt, copy the Systat13ENU32.msi (along with SYSTAT~1.cab, Qualit~1.cab, Monte_~1.cab and ExactT~1.cab) available in CD to C drive and run the below commands at the command prompt. (Please do not copy SilentInstallProperty.ini in the location).

```
msiexec /i C:\Systat13ENU32.msi SSI_SERIALNUMBER=xxxxxx
SSI_INSTTYPE=Network INSTALLDIR="D:\SYSTAT 13\" ALLUSERS=""
COMPANYNAME="My Company Name" USERNAME="My User Name" /qb
```

SYSTAT will be installed to all users by default. To install SYSTAT only for the current user on the same computer on NT/2K/XP operating systems use the property ALLUSERS="" in the command line as follows: This option can be combined with any of the above options.

```
msiexec /i C:\Systat13ENU32.msi SSI_SERIALNUMBER=xxxxxx
ALLUSERS="" /qb
```

Please note from above command that SSI_SERIALNUMBER property is mandatory and other property values will be assigned to default values by installer. If property is not specified at the command prompt, installer will initialize default value. Default property values are as follows .

```
SSI_INSTTYPE=standalone
INSTALLDIR="C:\Program Files\SYSTAT\SYSTAT 13"
ALLUSERS=1
COMPANYNAME= Organization name registered to Operating System
USERNAME= User Name registered to Operating System
```

The user's name, company name and the serial number (printed at the back on the product CD's jewel case) can be provided as follows during silent installation. This information will be displayed in the product's Help=>About box if provided.

```
msiexec /i "C:\Systat13ENU32.msi" USERNAME="John Doe"
COMPANYNAME="Your company name"
SSI_SERIALNUMBER="XXXXXXXXXX" /qb
```

Note: For silently installing Systat13 64 bit, use Systat13ENU64.msi and data.cab files by following above procedures.

Note: To install Systat13 application in Windows Vista and Windows7 operating systems silently, open the command prompt and change the directory to Systat13ENU32 folder CD (assume that CD ROM belongs F drive) and execute below command

```
F:\Systat13ENU32.exe /s /v/qb
```

Uninstalling SYSTAT

SYSTAT 13 can be uninstalled in any Windows environment using either of the following two methods:

1. Click *Start => Programs => SYSTAT 13 => Uninstall* to uninstall the program.
2. From the Windows **Start** menu choose:
 - **Settings => Control Panel**
 - Select **Add/Remove Programs**.
 - Choose **SYSTAT 13** from the list
 - Click **(Change/) Remove**.
 - Follow the instructions on the screen to remove SYSTAT.

To remove the product silently use the following command

```
msiexec /x C:\Systat13ENU32\Systat13ENU32.msi
```

Note: If SYSTAT is installed for all users then uninstalling SYSTAT for any user will uninstall SYSTAT for all users.

System Requirements

The minimum hardware and software requirements for SYSTAT are:

- Windows 32 bit (XP, Vista, 2003 Server and Windows 7, Windows 8.1, Windows 10)
- Windows 64 bit (XP, Vista, 2003 Server and Windows 7, Windows 8.1, Windows 10)
- Pentium level processor
- 512 MB RAM (1GB recommended)
- 600 MB Free Disk Space (inclusive of Manual PDFs ~87 MB and MSI Cache)
- CD-ROM Drive
- SVGA Adapter and Monitor
- Internet Explorer 7 (Internet Explorer 8 recommended)